



Determination of lithium in human serum by isotope dilution atomic absorption spectrometry

Alexander Winkelmann^{1,2} · Dalia Morcillo^{1,2} · Silke Richter² · Sebastian Recknagel² · Jens Riedel² · Jochen Vogl² · Ulrich Panne^{1,2} · Carlos Abad²

Received: 1 March 2021 / Revised: 15 July 2021 / Accepted: 26 August 2021 / Published online: 10 September 2021
© The Author(s) 2021

Abstract

The therapeutic dose of lithium (Li) compounds, which are widely used for the treatment of psychiatric and hematologic disorders, is close to its toxic level; therefore, drug monitoring protocols are mandatory. Herein, we propose a fast, simple, and low-cost analytical procedure for the traceable determination of Li concentration in human serum, based on the monitoring of the Li isotope dilution through the partially resolved isotope shift in its electronic transition around 670.80 nm using a commercially available high-resolution continuum source graphite furnace atomic absorption spectrometer. With this technique, serum samples only require acidic digestion before analysis. The procedure requires three measurements—an enriched ⁶Li spike, a mixture of a certified standard solution and spike, and a mixture of the sample and spike with a nominal ⁷Li/⁶Li ratio of 0.82. Lanthanum has been used as an internal spectral standard for wavelength correction. The spectra are described as the linear superposition of the contributions of the respective isotopes, each consisting of a spin-orbit doublet, which can be expressed as Gaussian components with constant spectral position and width and different relative intensity, reflecting the isotope ratio in the sample. Both the spectral constants and the correlation between isotope ratio and relative band intensity have been experimentally obtained using commercially available materials enriched with Li isotopes. The Li characteristic mass (m_c) obtained corresponds to 0.6 pg. The procedure has been validated using five human serum certified reference materials. The results are metrologically comparable and compatible to the certified values. The measurement uncertainties are comparable to those obtained by the more complex and expensive technique, isotope dilution mass spectrometry.

Keywords Lithium · Human serum · Isotope dilution · Atomic absorption spectrometry · High-resolution continuum source graphite furnace atomic absorption spectrometry

Introduction

The therapeutic effects of lithium (Li) on the brain and blood are well known and applied for the treatment of manic depression, granulocytopenia generated by radiation and chemotherapy, and immunoglobulin stimulation. Li has a narrow therapeutic range of 0.6–1.0 mmol L⁻¹. Concentrations of Li above 1.5 mmol L⁻¹ in serum can be toxic [1]; hence, close

monitoring of the Li concentration is mandatory to ensure effective and secure treatment. Considering the interindividual variation of the half-life of Li in patients, a quick and accurate analysis is required [2]. Currently, the primary analytical techniques used for Li quantification are based on relative measurements, which require reference serum standards for calibration, thereby hampering traceability and comparability.

Several analytical methods for the quantification of Li in biological samples have been proposed, including atomic absorption spectrophotometry (AAS) [3], inductively coupled plasma (ICP) optical emission spectroscopy (OES) [4], and ICP mass spectrometry (MS) [5]. The development of ICP-based instrumentation has facilitated the simultaneous quantification of several elements of medical interest with sufficient accuracy (trueness and precision) [6–8]. However, ICP-based methods generally require laborious sample preparation, and are expensive and matrix dependent. In this

Published in the topical collection celebrating *ABCs 20th Anniversary*.

✉ Carlos Abad
carlos.abad@bam.de

¹ Department of Chemistry, Humboldt-Universität zu Berlin, Brook-Taylor-Str. 2, 12489 Berlin, Germany

² Bundesanstalt für Materialforschung und -prüfung (BAM), Richard-Willstätter-Str. 11, 12489 Berlin, Germany

context, methodologies like isotope dilution (ID) are less matrix sensitive and provide reduced uncertainties in the measurement [9, 10]. ID is an absolute technique that improves metrological quality and provides an almost standard-free calibration approach, once the isotope spike is characterized.

Quantification of Li by ID-ICP-MS has been proposed for geological and biological samples [11, 12]. ID-MS can provide analyte concentrations that are traceable to the international system of units (SI) and thus enabling comparability, when certain requirements are fulfilled [13]. However, in the case of Li, MS-based methods are prone to large mass bias and matrix effects and huge instrumental drift, which is primarily due to a low mass-to-charge ratio and the resulting high relative mass difference between the two Li isotopes [14]. The high extent of the mass bias and its high fluctuation hinder its accurate correction, which in turn negatively affects the measurement uncertainty.

Some attractive alternatives to ID-MS techniques are based on the isotope shift in atomic and molecular electronic spectra. For example, ID-AAS was proposed by Brost et al. to monitor the absorption coefficient of naturally occurring Li in human plasma and enriched isotopic materials using natural and ^6Li -enriched hollow cathode lamps. The bias for Li recovery in plasma ranged from -2.8 to 0.6% [15]. However, the availability of isotopically enriched hollow cathode lamps limits the practical application of this method.

In recent years, optical spectrometry has flourished for isotope analysis. The use of a continuum light source coupled with a high-resolution echelle spectrometer and a charge-coupled device detector allows the monitoring of the electronic transition of transient diatomic molecules. The isotope shift in the electronic spectra of these diatomic molecules depends on the reduced mass and is larger than those observed in atomic spectra [16]. Thus, the monitoring of the isotope shift and the relative intensities of the isotopologue couple $\text{Al}^{35}\text{Cl}/\text{Al}^{37}\text{Cl}$ proved useful for the trace analysis of Cl via ID using a commercially available high-resolution continuum source graphite furnace molecular absorption spectrometer (HR-CS-GF-MAS) [17]. The same approach was applied for the trace analysis of Br and Ca [18, 19]. Additionally, accurate isotope ratios, which are the core of ID-MS applications, can be achieved by HR-CS-MAS and HR-CS-AAS using modern methods for data analysis like partial least square regression and machine learning, as was demonstrated for B and Li, respectively [20, 21]. In this latter work on Li isotope ratio analysis, the electronic transition $2^2\text{P} \leftarrow 2^2\text{S}$ was explored. This transition can be used for Li quantification by ID.

This work investigates the determination of the amount of Li in human serum reference materials based on an ID approach applied to HR-CS-GFAAS. Since the atomic spectra of Li follow a nonlinear Gaussian function, this can be used for isotope ratio analysis for SI traceability instead of a calibration model by machine learning. The performance

achieved by ID-HR-GFAAS in the present study is compared with those obtained by certificates of analysis of reference materials and previous reports on ID-MS.

Materials and methods

Sample preparation

High-purity deionized water with a resistivity of $18\text{ M}\Omega\text{ cm}$ obtained from a Milli-Q system (Millipore gradient, Merck Millipore, Darmstadt, Germany) was used throughout the experiments. Nitric acid (HNO_3 ; EMSURE®, Merck, Darmstadt, Germany) was used after purification by subboiling distillation in PFA containers. The samples analyzed consisted of five serum certified reference materials (CRMs): BCR 304 (Joint Research Centre, Belgium), ERM-DA250a and ERM-DA251a (LGC Limited, UK), Seronorm L-1 (SERO AS, Norway), and Seronorm L-2 (SERO AS, Norway). A spike solution was prepared from metallic ^6Li ($\approx 95\%$ ^6Li) in 2% HNO_3 with a mass concentration of around 0.5 mg L^{-1} . For the ID analysis, 0.3 g of each serum sample was digested in triplicate using 2 mL HNO_3 (35%) and 1 mL 15% H_2O_2 (Suprapur®, Merck, Darmstadt, Germany) at $100\text{ }^\circ\text{C}$ for 30 min . After digestion, 0.7 g of spike solution was added, and the solution was filled up to 25 mL with 2% HNO_3 . For reverse ID [9], 0.3 g of an ICP standard solution (traceable to NIST SRM® 3129a, Certipur®, Merck, Darmstadt, Germany) was mixed with 0.7 g of the spike solution and filled up to 10 mL with 2% HNO_3 . An ^7Li -enriched solution was used for optimization of the fit parameters and it was prepared by diluting $^7\text{Li}_2\text{CO}_3$ ($\geq 99\text{ atom}\%$, 99% ^7Li , Sigma-Aldrich, USA) in 2% HNO_3 for a mass concentration of around 0.5 mg L^{-1} .

Atomic absorption spectrometry measurements

A ContrAA 800D HR-CS-GFAAS model (Analytik Jena, Germany) with a graphite furnace (PIN platform) was used for all measurements. The wavelength of the instrument was centered at 670.7845 nm . The optimized measurement conditions provided in our recent work were applied [21]. However, owing to the low vapor pressure of La, the atomization temperature was increased to $2500\text{ }^\circ\text{C}$. For each measurement, $10\text{ }\mu\text{L}$ of digested sample was injected. The Li concentration of the injected solutions was adjusted to match an extinction of 3.0 integrated from 670.7361 to 670.8680 nm . A La ICP standard solution (Certipur®, Merck, Darmstadt, Germany) was added as an internal spectral standard with a final mass concentration of 2 g L^{-1} . Each sample was measured ten times, and 150 spectra were recorded during each atomization. The same procedure was applied for the ^6Li spike solution and the spiked standard solution (used for reverse ID) and for the

⁷Li solution. Between each sample, a blank sample was measured five times.

Data analysis

For each measurement, 150 spectra were collected during the Li atomization and averaged. Spectral data were preprocessed and reduced using the MATLAB software (R2020a, The MathWorks Inc., USA). The spectral data of each measurement were compiled, transformed, and integrated from three-dimensional to two-dimensional spectra by converting them to their average. Finally, the area of the Li line was normalized to the unity. Using the *fminsearch* function in MATLAB, the spectra were fit to Eq. 1, a sum of four Gaussian functions, where *c_i* is the central wavelength, *w_i* is the spectral width, and *A_i* is the peak area.

$$f(x) = \sum_{i=1}^4 \frac{A_i}{w_i \sqrt{\pi/2}} \cdot e^{-2\left(\frac{x-c_i}{w_i}\right)^2} \quad (1)$$

The four individual Gaussian functions correspond to two spin-orbit split isotope transitions. This results in tight constraints of the 12 variables in Eq. 1, enabling a robust fitting. All four central wavelength positions *c_i* are only governed by the term energies and can be treated as constant. The same applies for the spectral widths *w_i*, which are primarily determined by the instrumental resolution. Since the four peaks reflect only two doublet transitions of the two isotopes, the ratio between two respective areas (e.g., *A₁* and *A₂*) is constant. Namely, the spectral displacement between *c₁* and *c₃* is the isotopic shift of 15.80 pm between ⁷Li and ⁶Li [21]. The distance between *c₁* and *c₂* or *c₃* and *c₄* is the spin-orbit splitting of 15.08 pm. For calibration of the absolute spectral position, the La signal was fit to a single Gaussian function. The La peak center (*c_{La}*) was then used to correct the Li peak center. Therefore, *A₁* and *A₃* are the only free variables. The scripts used for data preprocessing and analysis are provided in the Supplementary information (ESM1 and ESM2).

Results

The Li characteristic mass (*m_c*) which represents 1% transmittance was found to be 0.6 pg. The ID analysis was based on the isotope shift of the electronic transition 2²P ← 2²S, which exhibits two spin-orbit components for each isotope, as we previously described [21]. The nonlinear fit parameters were optimized using the spectra of ⁷Li- and ⁶Li-enriched materials. The wavelength of the spectra was corrected by using La as an internal spectral standard for the electronic transition ⁴F_{5/2} ← ⁴G^o_{7/2} as shown in Fig. 1 [22].

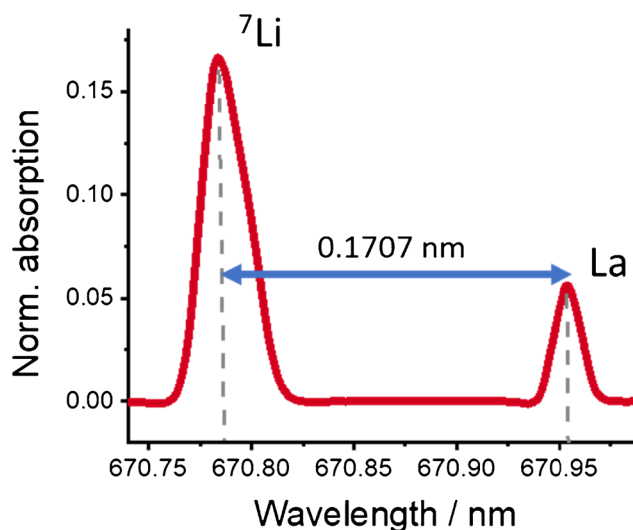


Fig. 1 Use of La as an internal spectral standard in a ⁷Li-enriched sample

The *A₁/A₂* ratio, which is dictated by spin statistics and the Einstein coefficients of the transitions, was determined to be 2.5645 (as well as the *A₃/A₄* ratio) and was fixed at this value. The widths *w₁* and *w₃* were determined and set to be 0.016 nm, and the widths *w₂* and *w₄* were fixed at 0.014 nm. For data analysis without correction, the peak center *c₁* was set at 670.7833 nm. The La-corrected Li peak center was set at *c₁* = *c_{La}* - 0.1707 nm. The ⁷Li/⁶Li isotope ratios were determined as the ratio of *A₁/A₃*. Subsequently, these ratios were fed into Eq. 2 [9] to calculate the Li contents in the serum samples, where *w_{sample}* is the Li mass fraction in the sample, *w_{spike,6Li}* is the ⁶Li mass fraction in the sample, *M(sample)* is the molar mass of Li in the sample, *M(⁶Li)* is the molar mass of ⁶Li, *x_{sample,6Li}* is the amount fraction of Li in the sample, *m* is the mass of sample and added spike, and *R* represents the isotope ratio in spike, sample, and spike-sample mix, respectively. This procedure for Li quantification via ID is illustrated in Fig. 2. Natural occurring variations in the Li isotope ratios are much smaller than the uncertainties in the ID-AAS measurements. Therefore, the atomic weight interval of Li with naturally isotopic composition, i.e., between 6.9387 and 6.9438 g mol⁻¹, was assumed in the serum samples [23]. The mass fraction of the spike solution was determined using reverse ID by rearranging the equation for *w_{spike}* and setting *w_{spike,6Li}* as the certified concentration of the ICP standard solution.

$$w_{sample} = w_{spike,6Li} \frac{M(sample)}{M(^6Li)} \cdot \frac{m_{spike}}{m_{sample}} \cdot \frac{R_{spike} - R_{mix}}{R_{mix} - R_{sample}} \quad (2)$$

Table 1 summarizes the results of the Li quantification by ID-HR-CS-GFAAS. The corresponding uncertainties were calculated according to the *Guide to the Expression of*

Table 1 Results of isotope dilution high-resolution continuum source graphite furnace atomic absorption spectrometry (ID-HR-CS-GFAAS) and certified Li content in human serum certified reference materials (CRMs)

CRM	Certified Li content		ID-AAS without correction			ID-AAS with La correction		
	Molarity (mmol L ⁻¹ Li)	<i>U</i>	Molarity (mmol L ⁻¹ Li)	<i>U</i>	<i>E_n</i>	Molarity (mmol L ⁻¹ Li)	<i>U</i>	<i>E_n</i>
BCR 304	0.985	0.029	0.92	0.09	0.69	0.976	0.022	0.25
ERM-DA250a	0.98*	0.06*	0.94	0.09	0.37	0.983	0.018	0.05
ERM-DA251a	0.66*	0.04*	0.62	0.07	0.50	0.663	0.022	0.07
Seronorm L1	0.73	0.15	0.82	0.06	0.57	0.830	0.016	0.68
Seronorm L2	1.45	0.29	1.52	0.13	0.22	1.463	0.028	0.04

The expanded uncertainties are expressed at the 95% confidence level ($k = 2$)

*Calculated from the certified mass fraction and its density. ERM-DA250a: $1.0294 \pm 0.0014 \text{ kg L}^{-1}$, ERM-DA251a: $1.0175 \pm 0.0013 \text{ kg L}^{-1}$. Atomic weight of lithium: $6.941 \pm 0.003 \text{ g mol}^{-1}$ [25]

Uncertainty in Measurement supplement 1 (GUM-S1) using *GUM Workbench 2.4* (Metrodata GmbH, Germany) and *Open Monte Carlo Engine v1.2.3* (Ruediger Kessel, NIST, USA). The metrological compatibility of the ID-HR-CS-GFAAS data with the certified values of the CRMs can be evaluated by applying the E_n value, which is the difference of two values divided by the expanded uncertainty of this difference [24]. Both values are considered metrologically compatible for $E_n < 1$ and not metrologically compatible for $E_n > 1$. As shown in Table 1, the Li concentrations determined by the present ID-HR-CS-GFAAS method are metrologically compatible with the certified values.

Discussion

A practical approach for Li quantification in human serum is presented, which is based on the partial resolution of the isotopic shift for the electronic transition $2^2\text{P} \leftarrow 2^2\text{S}$ around the wavelength of 670.80 nm by a commercially available HR-CS-GFAAS instrument. Monitoring of this isotope shift using the La peak as a reference enables the quantification of the isotopic components with high reproducibility. Therefore,

it can be used for an ID analysis by nonlinear fitting of four Gaussian curves representing the spin-orbit split doublets $^6\text{Li}_{1/2}$, $^6\text{Li}_{3/2}$, $^7\text{Li}_{1/2}$, and $^7\text{Li}_{3/2}$, respectively. The characterization of the individual transition properties (a doublet for each isotope) can be easily performed using commercially available ^6Li - and ^7Li -enriched materials. Although the HR-CS-GFAAS instrument provides high reproducibility for the wavelength measurements (± 0.5 pixels for resolution $\approx 140,000 \Delta\lambda/\lambda$), the accuracy is further improved by adding La as an internal spectral standard for wavelength correction. After characterization of the Gaussian parameters, only three samples need to be measured for Li quantification: (i) a ^6Li spike, (ii) a mixture of the standard and ^6Li spike, and (iii) a mixture of the sample and ^6Li spike.

Additionally, the human serum samples only require acidic dilution. No matrix effects are observed, at least for the type of biological samples evaluated. The application of this procedure for Li quantification in other relevant samples having complex matrices, such as geological samples and Li-ion batteries, is under investigation. The uncertainties (U) obtained using the present ID-HR-CS-GFAAS method with La correction lie between 0.016 and 0.028 mmol L⁻¹. The results are metrologically compatible with the certified reference material

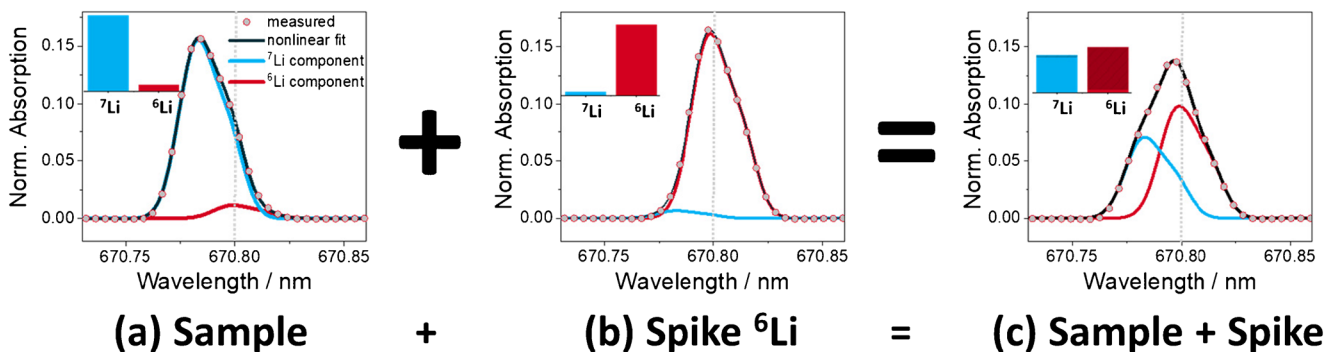


Fig. 2 Representation of the isotope dilution high-resolution continuum source graphite furnace atomic absorption spectrometry-based approach for determination of Li by its isotope shift in the atomic spectra. (a) Spectrum of a sample with a naturally occurring isotopic composition

(Li inductively coupled plasma standard) added to a ^6Li -enriched spike (b). The combination of (a) and (b) produced the spectrum of the mixture (c). On the top left corner of each spectrum, schematic representation of the isotope amount fractions of Li is provided

and comparable with those obtained with ID-AAS and ID-MS procedures, which afford uncertainties ranging between 0.006 and 0.015 mmol L⁻¹ [15, 26, 27]. The higher-order, traceability, and commutability of the studied certified reference materials assure their representability as clinical human serum samples [28–30]. For all these, ID-HR-CS-GFAAS can be envisaged as a metrological procedure for a fast, simple, and low-cost determination of Li in human serum.

Supplementary Information The online version contains supplementary material available at <https://doi.org/10.1007/s00216-021-03636-6>.

Acknowledgements The authors thank the support of Analytik Jena GmbH on this research with the HR-CS-AAS instrumentation.

Availability of data and material Data available on request from the authors.

Code availability The spectral data was analyzed using the MATLAB software (R2020a, The MathWorks Inc., USA). The scripts are provided in the Supplementary information.

Funding Open Access funding enabled and organized by Projekt DEAL.

Declarations

Ethics approval Not applicable.

Source of biological materials Human serum certified reference materials were acquired from commercial providers.

Consent to participate Not applicable.

Consent for publication The authors have approved the manuscript and agree with submission.

Conflict of interest The authors declare no competing interests.

Open Access This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons licence, and indicate if changes were made. The images or other third party material in this article are included in the article's Creative Commons licence, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons licence and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this licence, visit <http://creativecommons.org/licenses/by/4.0/>.

References

- Young W. Review of lithium effects on brain and blood. *Cell Transplant*. 2009;18(9):951–75. <https://doi.org/10.3727/096368909X471251>.
- Amdisen A. Lithium neurotoxicity—the reliability of serum lithium measurements. *Human Psychopharmacology: Clinical and Experimental*. 1990;5(3):281–5. <https://doi.org/10.1002/hup.470050324>.
- Manfro ID, Tegner M, Krutzmann ME, Artmann AC, Brandeburski MR, Peteffi GP, Linden R, Antunes MV. Determination of lithium in dried blood spots and dried plasma spots by graphite furnace atomic absorption spectrometry: method development, validation and clinical application. *Talanta*. 2020;216:120907. <https://doi.org/10.1016/j.talanta.2020.120907>.
- Iguchi K, Usuda K, Kono K, Dote T, Nishiura H, Shimahara M, Tanaka Y. Urinary lithium: distribution shape, reference values, and evaluation of exposure by inductively coupled plasma argon emission spectrometry. *J Anal Toxicol*. 1999;23(1):17–23. <https://doi.org/10.1093/jat/23.1.17>.
- Gouille JP, Mahieu L, Castermant J, Neveu N, Bonneau L, Laine G, Bouige D, Lacroix C. Metal and metalloid multi-elementary ICP-MS validation in whole blood, plasma, urine and hair. Reference values. *Forensic Sci Int*. 2005;153(1):39–44. <https://doi.org/10.1016/j.forsciint.2005.04.020>.
- Laur N, Kinscherf R, Pomytkin K, Kaiser L, Knes O, Deigner H-P. ICP-MS trace element analysis in serum and whole blood. *PLoS One*. 2020;15(5):e0233357. <https://doi.org/10.1371/journal.pone.0233357>.
- Forrer R, Gautschi K, Lutz H. Simultaneous measurement of the trace elements Al, As, B, Be, Cd, Co, Cu, Fe, Li, Mn, Mo, Ni, Rb, Se, Sr, and Zn in human serum and their reference ranges by ICP-MS. *Biol Trace Elem Res*. 2001;80(1):77. <https://doi.org/10.1385/BTER:80:1:77>.
- Bianchi F, Maffini M, Mangia A, Marengo E, Mucchino C. Experimental design optimization for the ICP-AES determination of li, Na, K, Al, Fe, Mn and Zn in human serum. *J Pharm Biomed Anal*. 2007;43(2):659–65. <https://doi.org/10.1016/j.jpba.2006.07.054>.
- Vogl J. Characterisation of reference materials by isotope dilution mass spectrometry. *J Anal At Spectrom*. 2007;22(5):475–92. <https://doi.org/10.1039/B614612K>.
- Brun V, Masselon C, Garin J, Dupuis A. Isotope dilution strategies for absolute quantitative proteomics. *J Proteome*. 2009;72(5):740–9. <https://doi.org/10.1016/j.jprot.2009.03.007>.
- Kramer U, Kress M, Reinauer H, Spannagl M, Kaiser P. Candidate reference measurement procedures for chloride, potassium, sodium, calcium, magnesium, and lithium by inductively coupled plasma (isotope dilution) sector field mass spectrometry (ICP-(ID) SFMS) in serum. *Clin Lab*. 2013;59(9–10):1017–29. <https://doi.org/10.7754/clin.lab.2012.120902>.
- Moriguti T, Makishima A, Nakamura E. Determination of Lithium contents in silicates by isotope dilution ICP-MS and its evaluation by isotope dilution thermal ionisation mass spectrometry. *Geostand Geoanal Res*. 2004;28(3):371–82. <https://doi.org/10.1111/j.1751-908X.2004.tb00756.x>.
- Milton MJT, Wielgosz RI. Uncertainty in SI-traceable measurements of amount of substance by isotope dilution mass spectrometry. *Metrologia*. 2000;37(3):199. <https://doi.org/10.1088/0026-1394/37/3/3>.
- Lin J, Liu Y, Hu Z, Yang L, Chen K, Chen H, Zong K, Gao S. Accurate determination of lithium isotope ratios by MC-ICP-MS without strict matrix-matching by using a novel washing method. *J Anal At Spectrom*. 2016;31(2):390–7. <https://doi.org/10.1039/c5ja00231a>.
- Brost DF, Brackett JM, Busch KW. Determination of lithium by optically monitored stable isotope dilution. *Anal Chem*. 1979;51(9):1512–6. <https://doi.org/10.1021/ac50045a035>.
- Russo RE, Bol'shakov AA, Mao X, McKay CP, Perry DL, Sorkhabi O. Laser ablation molecular isotopic spectrometry.

- Spectrochim Acta, Part B. 2011;66(2):99–104. <https://doi.org/10.1016/j.sab.2011.01.007>.
17. Nakadi FV, da Veiga MAMS, Aramendía M, García-Ruiz E, Resano M. Chlorine isotope determination via the monitoring of the AlCl molecule by high-resolution continuum source graphite furnace molecular absorption spectrometry – a case study. *J Anal At Spectrom.* 2015;30(7):1531–40. <https://doi.org/10.1039/c5ja00055f>.
 18. Nakadi FV, da Veiga MAMS, Aramendia M, Garcia-Ruiz E, Resano M. Br isotope determination via the monitoring of CaBr transitions using high-resolution continuum source graphite furnace molecular absorption spectrometry. Potential for direct determination of Br in solid samples using isotope dilution. *J Anal At Spectrom.* 2016;31(7):1381–90. <https://doi.org/10.1039/C6JA00114A>.
 19. Zanatta MBT, Nakadi FV, Resano M, da Veiga MAMS. Calcium isotope determination in urine samples via the monitoring of ^{44}CaF and ^{40}CaF molecules by high-resolution continuum source graphite furnace molecular absorption spectrometry. *J Anal At Spectrom.* 2019. <https://doi.org/10.1039/c9ja00233b>.
 20. Abad C, Florek S, Becker-Ross H, Huang M-D, Heinrich H-J, Recknagel S, Vogl J, Jakubowski N, Panne U. Determination of boron isotope ratios by high-resolution continuum source molecular absorption spectrometry using graphite furnace vaporizers. *Spectrochim Acta, Part B.* 2017;136(Supplement C):116–22. <https://doi.org/10.1016/j.sab.2017.08.012>.
 21. Winckelmann A, Nowak S, Richter S, Recknagel S, Riedel J, Vogl J, Panne U, Abad C. High-resolution atomic absorption spectrometry combined with machine learning data processing for isotope amount ratio analysis of lithium. *Anal Chem.* 2021. <https://doi.org/10.1021/acs.analchem.1c00206>.
 22. Hartog EAD, Palmer AJ, Lawler JE. Radiative lifetimes and transition probabilities of neutral lanthanum. *J Phys B Atomic Mol Phys.* 2015;48(15). <https://doi.org/10.1088/0953-4075/48/15/155001>.
 23. Laeter JR, Böhlke JK, Bièvre PD, Hidaka H, Peiser HS, Rosman KJR, Taylor PDP. Atomic weights of the elements. Review 2000 (IUPAC Technical Report). *Pure Appl Chem.* 2003;75(6):683–800. <https://doi.org/10.1351/pac200375060683>.
 24. Vogl J, Rosner M, Kasemann SA, Kraft R, Meixner A, Noordmann J, Rabb S, Rienitz O, Schuessler JA, Tatzel M, Vocke RD. Intercalibration of Mg isotope delta scales and realisation of Si traceability for Mg isotope amount ratios and isotope delta values. *Geostand Geoanal Res.* 2020;44(3):439–57. <https://doi.org/10.1111/ggr.12327>.
 25. Meija J, Coplen TB, Berglund M, Brand WA, De Bièvre P, Gröning M, Holden NE, Irrgeher J, Loss RD, Walczyk T, Prohaska T. Atomic weights of the elements 2013 (IUPAC Technical Report). *Pure Appl Chem.* 2016;88(3):265–91. <https://doi.org/10.1515/pac-2015-0305>.
 26. Michiels E, De Bièvre P. Accurate assay of lithium by isotope dilution mass spectrometry. *Int J Mass Spectrom.* 1983;48:369–72. [https://doi.org/10.1016/0020-7381\(83\)87104-6](https://doi.org/10.1016/0020-7381(83)87104-6).
 27. Lamberty A, De Bièvre P. CBNM - determination of Li in BCR reference materials RM 303 and RM 304 lyophilized serum by isotope dilution mass spectrometry (IDMS). *Fresenius J Anal Chem.* 1993;345(2):330–4. <https://doi.org/10.1007/BF00322624>.
 28. Infusino I, Frusciant E, Ferrero CA, Panteghini M. Commutability of two JCTLM-listed secondary reference materials for two commercial lithium assays. *Clin Chim Acta.* 2012;414:152–3. <https://doi.org/10.1016/j.cca.2012.09.001>.
 29. Barbagallo RP, Boley N, Holcombe G, Merson S, Mussell C, Pritchard C, Stokes P, Wood S, Ducroq D, Thomas A. Production and certification of four frozen human serum certified reference materials containing creatinine and electrolytes. *Ann Clin Biochem.* 2008;45(2):160–6. <https://doi.org/10.1258/acb.2007.007126>.
 30. Panteghini M, Braga F. Implementation of metrological traceability in laboratory medicine: where we are and what is missing. *Clinical Chemistry and Laboratory Medicine (CCLM).* 2020;58(8):1200–4. <https://doi.org/10.1515/cclm-2019-1128>.
- Publisher's note** Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.